

Dental Implant Unit

MOTORSURG

Always Reliable



SPECIFICATIONS:

Speed(rpm)	200-40,000
Torque(N·cm)	5-80
Max. pump output(ml/min)	130
Program	8
Screen	LED color touch
Foot pedal	IPX6
One-button calibration	Yes
Light	White LED
Length of motor cable(m)	2
Volume(cm)	24x26.5x15
Weight(kg)	3.65
Power supply	220-240V & 100-120V 50/60Hz




<https://stomshop.pro/>



POWERFUL & PRECISE

- German Bearing & Powerful Motor
5N.cm-80N.cm torque and 200rpm-40,000rpm speed provide high performance.

- Real-time Torque

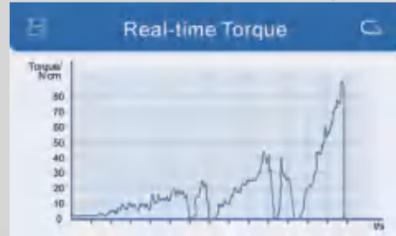

- One-button Calibration



- LED Contra Angle



SAFE & STABLE

- Stable under high torque
Ultra low speed available
- Surgical parameters can be preset to avoid misoperation.
- Smart error reminder, effectively guiding users to solve common problems in surgery.
- Visualization of surgical procedures, guide users to do surgery and improve the success rate.



Marking Pilot drilling Template drilling Tapping Placing implant Setting a closure cap Free use Rinsing



STRONG CONTRA ANGLE

- Adjustable brightness helps surgery more clearly and precisely.
- Put the water inlet on the bottom.




- * Safe: avoid scratches when disassembling.
- * Comfortable: there is no water pipe at the front of the contra angle, which takes up less space in the mouth and makes the patient more comfortable.
- Ergonomically designed and well-balanced
- Comfortable to hold without fatigue
- Detachable
- Durable




USER-FRIENDLY

- LED color touch screen: simple interface, intuitive, sensitive and easy to use

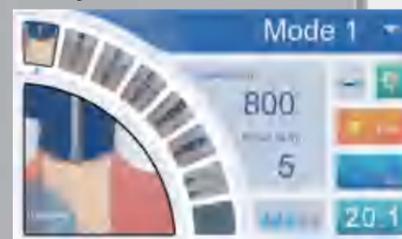

- Data recording
Data can be downloaded and analyzed by USB flash disk.

- 8 modes: enter the surgery quickly and dentist can set parameters and customize the mode.

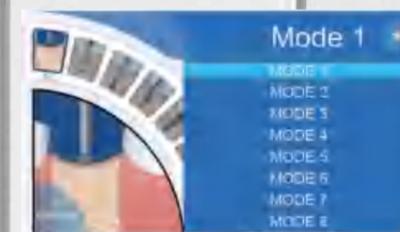

*USB flash disk isn't included.

MULTI-FUNCTIONAL FOOT PEDAL

- Intelligent, user-friendly and aseptic operation
- Compared with other brands, ours can be started by lightly stepping on it, which is more convenient.



- ① Start/stop button
- ② Water volume control button
- ③ Program switch button
- ④ Forward/reverse switch button



QUIET WATER PUMP

- High performance: maximum water volume can reach 130ml/min.
- Adjustable water volume
- Simple connection

