

SAM 3

The SAM 3 Articulator has all of the basic features of the SAM 2P and those of the SAM 2PX. Additionally, it has a unique and different high precision centric locking mechanism which functions independent

of both Bennett Guidance and Condylar Pathway Inserts. An optional micro adjustable protrusion and retrusion screw is available.



Image of
SAM 3 Professional

SAM 3 Professional

includes all 3 curvature inserts (1-3), all 4 Bennett guidance inserts (straight, small, extended and immediate sideshifts), adjustable incisal table, micro adjustable incisal pin, micro adjustable protrusion (to 6 mm) and retrusion (to 2 mm) screw, vertical support rods and tilt support rods

SAM 3 Standard

includes green curvature 2 insert, Bennett guidance insert (straight, metal), flat incisal table, micro adjustable protrusion (to 6 mm) and retrusion (to 2 mm) screw, vertical support rods and tilt support rods

SAM 3 Basic

includes green curvature 2 insert, Bennett guidance insert (straight, metal), flat incisal table


SAM 3 Features



Centric lock with cone ring



Bennett angulation indicator and setting mechanism

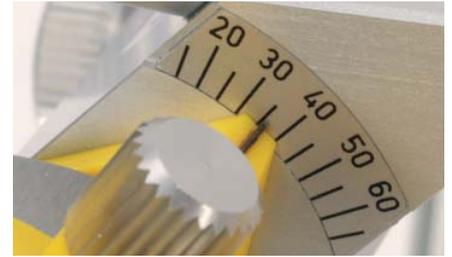


Condylar Housing inclination angle indicator



Micro Adjustable Incisal Pin (included in SAM 3 Professional)



Facebow earpiece pin

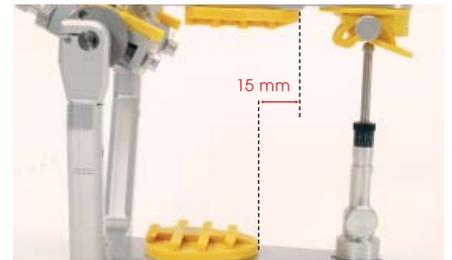

Mandibular mounting plate 15 mm sagittal posterior position



Selection of cast mounting plates, MPS magnetic or screw type. (AXIOSPLIT® optional)



Vertical Support Rod for upper member of articulator.



Tilt support rods for angular positioning of the articulator



Micro Adjustable Protrusion (to 6 mm) and Retrusion (to 2 mm) Screw (not included in SAM 3 Basic)



Transparent condylar housing cover



Adjustable Incisal Table (included in SAM 3 Professional)



Options: Bennett and Condylar Insert Set 3



Visibly high quality and precision assures reliability and longevity